

PRODUCT INFORMATION AND MANUAL

Mouse TNF- α ELISA

BMS607/2 BMS607/2TWO BMS607/2TEN

Enzyme-linked immunosorbent assay for quantitative

detection of mouse TNF- α .

For research use only.

Not for diagnostic or therapeutic procedures.

**Mouse TNF- α
BMS607/2
BMS607/2TWO
BMS607/2TEN**



Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe

TABLE OF CONTENTS

1	Intended Use	3
2	Summary	3
3	Principles of the Test	5
4	Reagents Provided	7
5	Storage Instructions – ELISA Kit	10
6	Specimen Collection and Storage Instructions	10
7	Materials Required But Not Provided	11
8	Precautions for Use	12
9	Preparation of Reagents	14
10	Test Protocol	19
11	Calculation of Results	24
12	Limitations	27
13	Performance Characteristics	28
14	Ordering Information	33
15	Reagent Preparation Summary	33
16	Test Protocol Summary	35

1 Intended Use

The mouse TNF- α ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of mouse TNF- α . **The mouse TNF- α ELISA is for research use only. Not for diagnostic or therapeutic procedures.**

2 Summary

TNF- α is a multifunctional cytokine involved in many different pathways, in homeostasis and pathophysiology of mammals. It can show opposing biological effects suggesting complex regulatory mechanisms.

TNF- α , also known as cachectin, was first detected as a cytotoxic factor inducing lysis of certain tumor cells. The TNF- α gene is member 2 of the TNF-superfamily (consisting of at least 20 distinct members).

TNF- α release is mainly triggered by viral infections and endotoxins, lipopolysaccharides or other bacterial components, by tissue injury, DNA-damage and by IL-1, PDGF and TNF- α itself. It is primarily expressed in macrophages, but also in monocytes, neutrophils, NK-cells, mast-cells, endothelial cells and activated lymphocytes. TNF- α expression in endothelial cells and fibroblasts can be induced by IL-17.

The expression of other cytokines, chemokines, reactive oxygen intermediates, nitric oxide and prostaglandins is stimulated by TNF- α . The initially membrane bound TNF- α is enzymatically cleaved by TACE (= ADAM17). The soluble monomers aggregate to homotrimers and are secreted into blood and other biological fluids. The membrane bound and the soluble form are biologically active and bind to the TNF-receptors TNFR1 (= TNFRSF1A, p55-60) and TNFR2 (= TNFRSF1B, TNFR2, p75-80).

Upon ligand binding, receptors form trimers leading to conformational changes, protein dissociation (SODD = silencer of death domains, BAG4, Bcl2-associated athanogene 4) and association (TRADD = TNF-R1 associated death domain protein) and yielding the following biological activities:

- transcription of anti-apoptotic factors and proteins involved in cell proliferation and inflammation via binding of TRAF2 (TNF-R associated factor 2) and RIPK1 (TNF-R interacting serine-threonine kinase 1) and activation of the transcription factor NF- κ B.
- cell proliferation, differentiation but also apoptosis via TRAF2 binding, kinases activation, activation of c-Jun and ATF2 (JNK-MAPK-pathway).

- apoptosis via the binding of FADD (Fas associating protein with death domain) to TRADD and activation of caspases (including caspase 8 = FLICE).

- necrosis, a caspase independent cell death, mediated by NADPH oxidases, which form a complex with TRADD and RIPK1, leading to the generation of oxygen species.

TNF-R2 contains no DD (death domain), but exhibits its function via direct TRAF binding.

Thus the multiple biological functions of TNF- α comprise cellular proliferation and differentiation, tumorigenesis, apoptotic or necrotic cell death (including certain tumor cell lines), immunoregulatory activities, lipid metabolism, coagulation and endothelial function. It promotes local or systemic inflammation (TNF- α is a potent pyrogen) and stimulates the acute phase response. Very high expressions of TNF- α after infection can lead to septic shock (TNF- α is highly cytotoxic.), whereas sustained low levels induce cachexia and inflammation.


Dysregulation of TNF- α is involved in many diseases.

For literature update refer to **www.bendermedsystems.com**

3 Principles of the Test


An anti-mouse TNF- α coating antibody is adsorbed onto microwells.

Figure 1

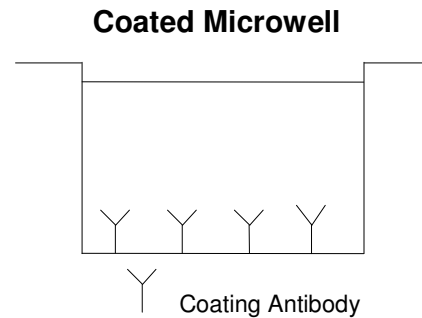


Mouse TNF- α present in the sample or standard binds to antibodies adsorbed to the microwells.

Figure 2

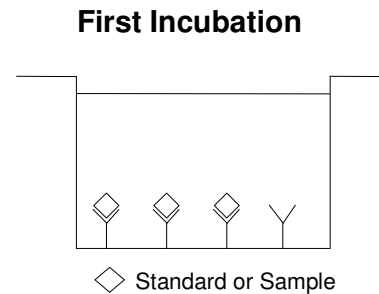


Following incubation unbound biological components are removed during a wash step and a biotin-conjugated anti-mouse TNF- α antibody is added and binds to mouse TNF- α captured by the first antibody.

Figure 3

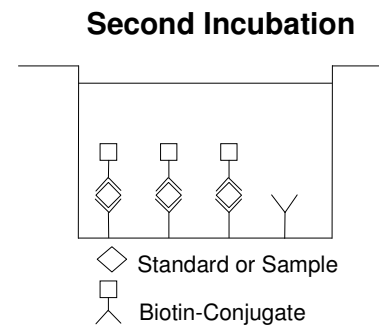

Following incubation unbound biotin-conjugated anti-mouse TNF- α antibody is removed during a wash step. Streptavidin-HRP is added and binds to the biotin-conjugated anti-mouse TNF- α antibody.

Figure 4

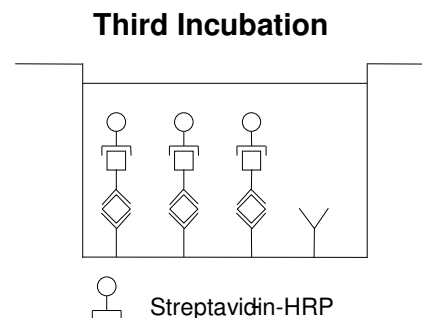

Following incubation unbound Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of mouse TNF- α present in the sample or standard. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from 7 mouse TNF- α standard dilutions and mouse TNF- α sample concentration determined.

Figure 5


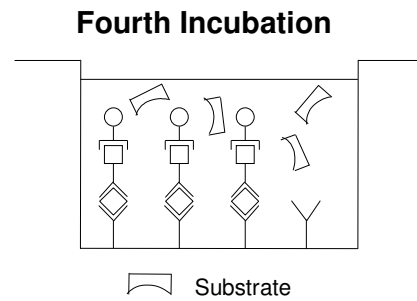


Figure 6

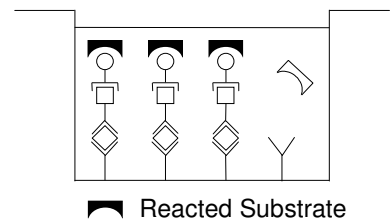

4 Reagents Provided

4.1 Reagents for mouse TNF- α ELISA BMS607/2 (96 tests)

- 1 aluminium pouch with a **Microwell Plate coated** with polyclonal antibody to mouse TNF- α
- 1 vial (150 μ l) **Biotin-Conjugate** anti-mouse TNF- α polyclonal antibody
- 1 vial (150 μ l) **Streptavidin-HRP**
- 2 vials mouse TNF- α **Standard** lyophilized, 2 ng/ml upon reconstitution
- 1 vial (12 ml) **Sample Diluent**
- 1 vial (5 ml) **Assay Buffer Concentrate 20x**
(PBS with 1% Tween 20 and 10% BSA)
- 1 bottle (50 ml) **Wash Buffer Concentrate 20x** (PBS with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 1 vial (0.4 ml) **Blue-Dye**
- 1 vial (0.4 ml) **Green-Dye**
- 1 vial (0.4 ml) **Red-Dye**
- 6 **Adhesive Films**

4.2 Reagents for mouse TNF- α ELISA BMS607/2TWO (2x96 tests)

- 2 aluminium pouches with a **Microwell Plate coated** with polyclonal antibody to mouse TNF- α
- 2 vials (150 μ l) **Biotin-Conjugate** anti-mouse TNF- α polyclonal antibody
- 2 vials (150 μ l) **Streptavidin-HRP**
- 4 vials mouse TNF- α **Standard** lyophilized, 2 ng/ml upon reconstitution
- 2 vials (12 ml) **Sample Diluent**
- 1 vial (5 ml) **Assay Buffer Concentrate 20x**
(PBS with 1% Tween 20 and 10% BSA)
- 2 bottles (50 ml) **Wash Buffer Concentrate 20x**
(PBS with 1% Tween 20)
- 2 vials (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 2 vials (15 ml) **Stop Solution** (1M Phosphoric acid)
- 1 vial (0.4 ml) **Blue-Dye**
- 1 vial (0.4 ml) **Green-Dye**
- 1 vial (0.4 ml) **Red-Dye**
- 12 **Adhesive Films**

4.3 Reagents for mouse TNF- α ELISA BMS607/2TEN (10x96 tests)

- 10 aluminium pouches with a **Microwell Plate coated** with polyclonal antibody to mouse TNF- α
- 10 vials (150 μ l) **Biotin-Conjugate** anti-mouse TNF- α polyclonal antibody
- 10 vials (150 μ l) **Streptavidin-HRP**
- 10 vials mouse TNF- α **Standard** lyophilized, 2 ng/ml upon reconstitution
- 10 vials (12 ml) **Sample Diluent**
- 2 vials (5 ml) **Assay Buffer Concentrate 20x**
(PBS with 1% Tween 20 and 10% BSA)
- 7 bottles (50 ml) **Wash Buffer Concentrate 20x**
(PBS with 1% Tween 20)
- 10 vials (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 10 vials (15 ml) **Stop Solution** (1M Phosphoric acid)
- 6 vials (0.4 ml) **Blue-Dye**
- 6 vials (0.4 ml) **Green-Dye**
- 6 vials (0.4 ml) **Red-Dye**
- 30 **Adhesive Films**

5 Storage Instructions – ELISA Kit

Store kit reagents between 2° and 8°C. Immediately after use remaining reagents should be returned to cold storage (2° to 8°C). Expiry of the kit and reagents is stated on labels.

Expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

6 Specimen Collection and Storage Instructions

Cell culture supernatant, serum and plasma (EDTA, citrate) were tested with this assay. Other biological samples might be suitable for use in the assay. Remove serum or plasma from the clot or cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20°C to avoid loss of bioactive mouse TNF- α . If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13.5).

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- 50 μ l to 300 μ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.

- Exposure to acid inactivates the conjugate.
- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents

Buffer Concentrates should be brought to room temperature and should be diluted before starting the test procedure. If crystals have formed in the **Buffer Concentrates**, warm them gently until they have completely dissolved.

9.1 Wash Buffer (1x)

Pour entire contents (50 ml) of the **Wash Buffer Concentrate** (20x) into a clean 1000 ml graduated cylinder. Bring to final volume of 1000 ml with glass-distilled or deionized water. Mix gently to avoid foaming. The pH of the final solution should adjust to 7.4.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer (1x) is stable for 30 days.

Wash Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Wash Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

9.2 Assay Buffer (1x)

Pour the entire contents (5 ml) of the **Assay Buffer Concentrate** (20x) into a clean 100 ml graduated cylinder. Bring to final volume of 100 ml with distilled water. Mix gently to avoid foaming.

Store at 2° to 8°C. Please note that the Assay Buffer (1x) is stable for 30 days.

Assay Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Assay Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

9.3 Biotin-Conjugate

Please note that the Biotin-Conjugate should be used within 30 minutes after dilution.

Make a 1:100 dilution of the concentrated **Biotin-Conjugate** solution with Assay Buffer (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	Biotin-Conjugate (ml)	Assay Buffer (1x) (ml)
1 - 6	0.06	5.94
1 - 12	0.12	11.88

9.4 Streptavidin-HRP

Please note that the Streptavidin-HRP should be used within 30 minutes after dilution.

Make a 1:200 dilution of the concentrated **Streptavidin-HRP** solution with Assay Buffer (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	Streptavidin-HRP (ml)	Assay Buffer (1x) (ml)
1 - 6	0.03	5.97
1 - 12	0.06	11.94

9.5 Mouse TNF- α Standard

Reconstitute **mouse TNF- α standard** by addition of distilled water. Reconstitution volume is stated on the label of the standard vial. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 2 ng/ml). Allow the reconstituted standard to sit for a 10-30 minutes. Mix well prior to making dilutions.

After usage remaining standard cannot be stored and has to be discarded.

Standard dilutions can be prepared directly on the microwell plate (see 10.c) or alternatively in tubes (see 9.5.1).

9.5.1 External Standard Dilution

Label 7 tubes, one for each standard point.

S1, S2, S3, S4, S5, S6, S7

Then prepare 1:2 serial dilutions for the standard curve as follows:

Pipette 225 μ l of Sample Diluent into each tube.


Pipette 225 μ l of reconstituted standard (concentration of standard = 2 ng/ml) into the first tube, labelled S1, and mix (concentration of standard 1 = 1 ng/ml).

Pipette 225 μ l of this dilution into the second tube, labelled S2, and mix thoroughly before the next transfer.

Repeat serial dilutions 5 more times thus creating the points of the standard curve (see Figure 7).

Sample Diluent serves as blank.

Figure 7

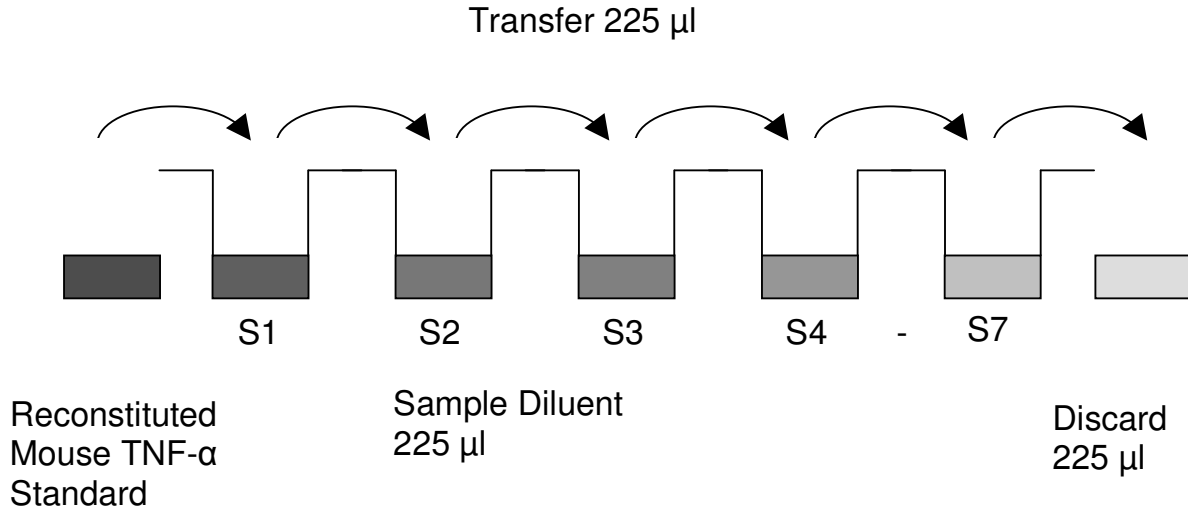

9.6 Addition of Colour-giving Reagents: Blue-Dye, Green-Dye, Red-Dye

In order to help our customers to avoid any mistakes in pipetting the Bender MedSystems ELISAs, Bender MedSystems offers a tool that helps to monitor the addition of even very small volumes of a solution to the reaction well by giving distinctive colours to each step of the ELISA procedure.

This procedure is optional, does not in any way interfere with the test results, and is designed to help the customer with the performance of the test, but can also be omitted, just following the instruction booklet.

Alternatively, the dye solutions from the stocks provided (**Blue-Dye, Green-Dye, Red-Dye**) can be added to the reagents according to the following guidelines:

1. Diluent:

Before standard dilution add the **Blue-Dye** at a dilution of 1:250 (see table below) to the appropriate diluent (1x) according to the test protocol. After addition of **Blue-Dye**, proceed according to the instruction booklet.

5 ml Sample Diluent	20 μ l Blue-Dye
12 ml Sample Diluent	48 μ l Blue-Dye
50 ml Sample Diluent	200 μ l Blue-Dye

2. Biotin-Conjugate: Before dilution of the concentrated Biotin-Conjugate, add the **Green-Dye** at a dilution of 1:100 (see table below) to the Assay Buffer (1x) used for the final conjugate dilution. Proceed after addition of **Green-Dye** according to the instruction booklet: Preparation of Biotin-Conjugate.

3 ml Assay Buffer (1x)	30 µl Green-Dye
6 ml Assay Buffer (1x)	60 µl Green-Dye
12 ml Assay Buffer (1x)	120 µl Green-Dye


3. Streptavidin-HRP: Before dilution of the concentrated Streptavidin-HRP, add the **Red-Dye** at a dilution of 1:250 (see table below) to the Assay Buffer (1x) used for the final Streptavidin-HRP dilution. Proceed after addition of **Red-Dye** according to the instruction booklet: Preparation of Streptavidin-HRP.

6 ml Assay Buffer (1x)	24 µl Red-Dye
12 ml Assay Buffer (1x)	48 µl Red-Dye

10 Test Protocol

- a. Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Each sample, standard, blank and optional control sample should be assayed in duplicate. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2°-8°C sealed tightly.
- b. Wash the microwell strips twice with approximately 400 µl **Wash Buffer** per well with thorough aspiration of microwell contents between washes. Allow the Wash Buffer to sit in the wells for about **10 – 15 seconds** before aspiration. Take care not to scratch the surface of the microwells.
After the last wash step, empty wells and tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing. Alternatively microwell strips can be placed upside down on a wet absorbent paper for not longer than 15 minutes. **Do not allow wells to dry.**
- c. **Standard dilution on the microwell plate** (Alternatively the standard dilution can be prepared in tubes - see 9.5.1):
Add 100 µl of Sample Diluent in duplicate to all **standard wells**. Pipette 100 µl of prepared **standard** (see Preparation of Standard 9.5, concentration = 2000.0 pg/ml) in duplicate into well A1 and A2 (see Table 1). Mix the contents of wells A1 and A2 by repeated aspiration and ejection (concentration of standard 1, S1 = 1000.0 pg/ml), and transfer 100 µl to wells B1 and B2, respectively (see Figure 8). Take care not to scratch the inner surface of the microwells. Continue this procedure 5 times, creating two rows of mouse TNF-α standard dilutions ranging from 1000.0 to 15.6 pg/ml. Discard 100 µl of the contents from the last microwells (G1, G2) used.

Figure 8

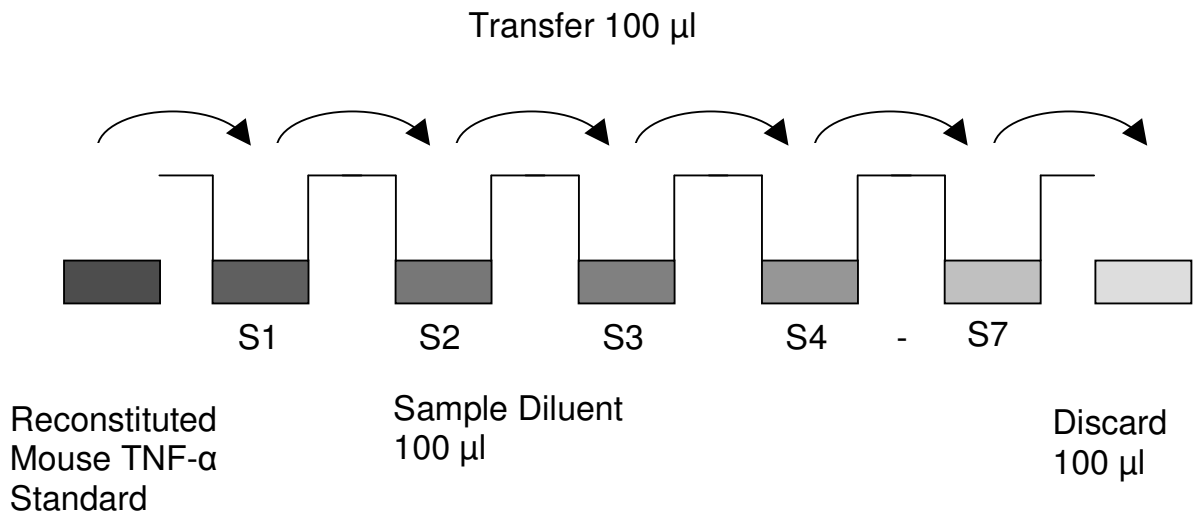

In case of an **external standard dilution** (see 9.5.1), pipette 100 μ l of these standard dilutions (S1 - S7) in the standard wells according to Table 1.

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (1000.0 pg/ml)	Standard 1 (1000.0 pg/ml)	Sample 1	Sample 1
B	Standard 2 (500.0 pg/ml)	Standard 2 (500.0 pg/ml)	Sample 2	Sample 2
C	Standard 3 (250.0 pg/ml)	Standard 3 (250.0 pg/ml)	Sample 3	Sample 3
D	Standard 4 (125.0 pg/ml)	Standard 4 (125.0 pg/ml)	Sample 4	Sample 4
E	Standard 5 (62.5 pg/ml)	Standard 5 (62.5 pg/ml)	Sample 5	Sample 5
F	Standard 6 (31.3 pg/ml)	Standard 6 (31.3 pg/ml)	Sample 6	Sample 6
G	Standard 7 (15.6 pg/ml)	Standard 7 (15.6 pg/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- d. Add 100 μ l of **Sample Diluent** in duplicate to the **blank wells**.
- e. Add 50 μ l of **Sample Diluent** to the **sample wells**.
- f. Add 50 μ l of each **sample** in duplicate to the **sample wells**.
- g. Cover with an adhesive film and incubate at room temperature (18 to 25°C) for 2 hours, if available on a microplate shaker set at 100 rpm.
- h. Prepare **Biotin-Conjugate** (see Preparation of Biotin-Conjugate 9.3).
- i. Remove adhesive film and empty wells. **Wash** microwell strips 5 times according to point b. of the test protocol. Proceed immediately to the next step.
- j. Add 100 μ l of **Biotin-Conjugate** to all wells.
- k. Cover with an adhesive film and incubate at room temperature (18 to 25°C) for 1 hour, if available on a microplate shaker set at 100 rpm.
- l. Prepare **Streptavidin-HRP** (refer to Preparation of Streptavidin-HRP 9.4).
- m. Remove adhesive film and empty wells. **Wash** microwell strips 5 times according to point b. of the test protocol. Proceed immediately to the next step.
- n. Add 100 μ l of diluted **Streptavidin-HRP** to all wells, including the blank wells.
- o. Cover with an adhesive film and incubate at room temperature (18° to 25°C) for 1 hour, if available on a microplate shaker set at 100 rpm.
- p. Remove adhesive film and empty wells. **Wash** microwell strips 5 times according to point b. of the test protocol. Proceed immediately to the next step.
- q. Pipette 100 μ l of **TMB Substrate Solution** to all wells.
- r. Incubate the microwell strips at room temperature (18° to 25°C) for about 10 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see next point of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. Alternatively the colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.9 – 0.95.

- s. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
- t. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean value.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the mouse TNF- α concentration on the abscissa. Draw a best fit curve through the points of the graph (a 5-parameter curve fit is recommended).
- To determine the concentration of circulating mouse TNF- α for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding mouse TNF- α concentration.
- **If instructions in this protocol have been followed samples have been diluted 1:2 (50 μ l sample + 50 μ l Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 2).**
- **Calculation of samples with a concentration exceeding standard 1 may result in incorrect, low mouse TNF- α levels. Such samples require further external predilution according to expected mouse TNF- α values with Sample Diluent in order to precisely quantitate the actual mouse TNF- α level.**
- It is suggested that each testing facility establishes a control sample of known mouse TNF- α concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 9. This curve cannot be used to derive test results. Each laboratory must prepare a standard curve for each group of microwell strips assayed.

Figure 9

Representative standard curve for mouse TNF- α ELISA. Mouse TNF- α was diluted in serial 2-fold steps in Sample Diluent. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

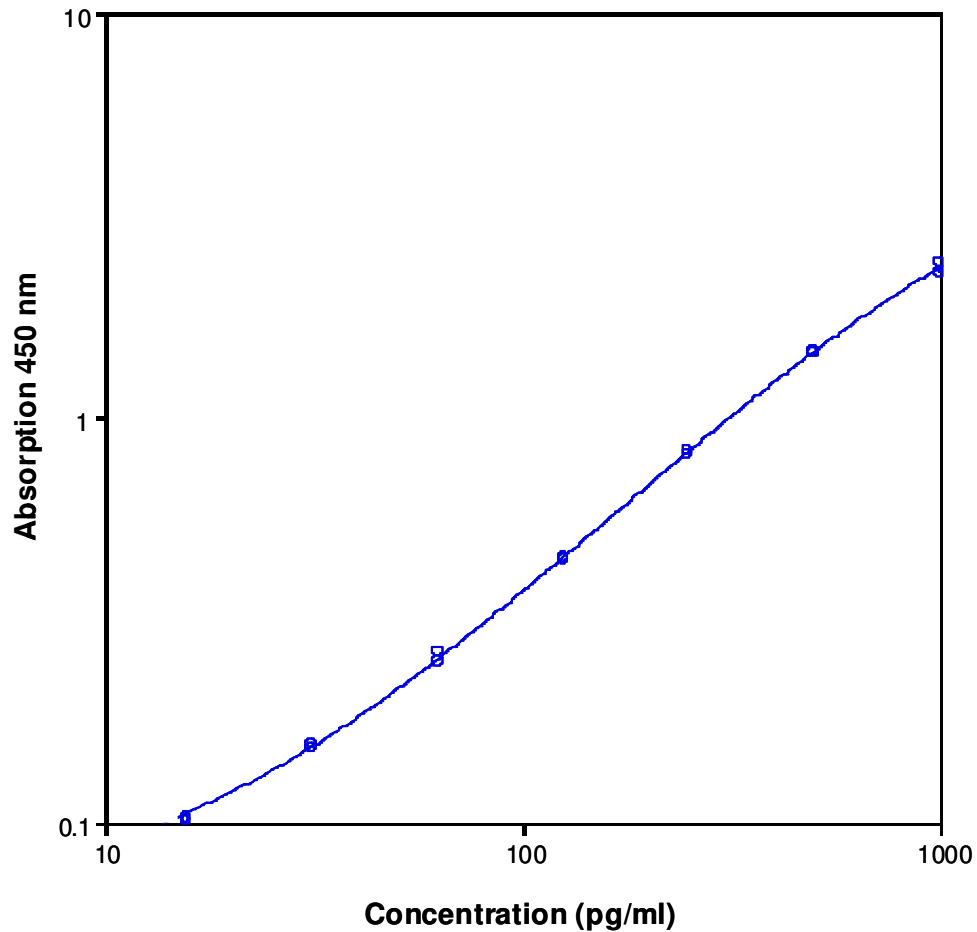


Table 2

Typical data using the mouse TNF- α ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	Mouse TNF- α Concentration (pg/ml)	O.D. at 450 nm	Mean O.D. at 450 nm	C.V. (%)
1	1000.0	2.296 2.403	2.349	2.3
2	500.0	1.449 1.467	1.458	0.6
3	250.0	0.830 0.815	0.823	0.9
4	125.0	0.452 0.447	0.450	0.5
5	62.5	0.265 0.251	0.258	2.6
6	31.3	0.157 0.154	0.156	1.1
7	15.6	0.103 0.101	0.102	1.2
Blank	0	0.040 0.041	0.040	

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of mouse TNF- α defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 3 pg/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 4 replicates of 8 serum samples containing different concentrations of mouse TNF- α . 2 standard curves were run on each plate. Data below show the mean mouse TNF- α concentration and the coefficient of variation for each sample (see Table 3). The calculated overall intra-assay coefficient of variation was 4.8%.

Table 3

The mean mouse TNF- α concentration and the coefficient of variation for each sample

Sample	Experiment	Mean Mouse TNF- α Concentration (pg/ml)	Coefficient of Variation (%)
1	1	422	6.4
	2	476	3.4
	3	432	2.8
2	1	1280	10.5
	2	1290	4.8
	3	1331	5.9
3	1	74	5.9
	2	73	6.2
	3	78	4.9
4	1	21	5.4
	2	20	8.4
	3	23	3.7
5	1	1271	4.2
	2	1389	3.5
	3	1348	5.0
6	1	474	2.2
	2	483	4.1
	3	458	5.0
7	1	262	3.0
	2	268	4.9
	3	253	4.6
8	1	217	3.6
	2	217	4.4
	3	194	2.7

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments. Each assay was carried out with 4 replicates of 8 serum samples containing different concentrations of mouse TNF- α . 2 standard curves were run on each plate. Data below show the mean mouse TNF- α concentration and the coefficient of variation calculated on 18 determinations of each sample (see Table 4). The calculated overall inter-assay coefficient of variation was 4.4%.

Table 4

The mean mouse TNF- α concentration and the coefficient of variation of each sample

Sample	Mean Mouse TNF- α Concentration (pg/ml)	Coefficient of Variation (%)
1	443	6.4
2	1300	2.1
3	75	3.3
4	22	7.2
5	1336	4.5
6	471	2.7
7	261	2.9
8	209	6.3

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of mouse TNF- α into serum. Recoveries were determined in 3 independent experiments with 4 replicates each. The unspiked serum was used as blank in these experiments.

The recovery ranged from 68% to 92% with an overall mean recovery of 80%.

13.4 Dilution Parallelism

Serum samples with different levels of mouse TNF- α were analysed at serial 2 fold dilutions with 4 replicates each.

The recovery ranged from 97.2% to 112.1% with an overall recovery of 104.0% (see Table 5).

Table 5

Sample	Dilution	Expected Mouse TNF- α Concentration (pg/ml)	Observed Mouse TNF- α Concentration (pg/ml)	Recovery of Expected Mouse TNF- α Concentration (%)
1	1:2		1422	
	1:4	711	712	100.1
	1:8	355	370	104.1
	1:16	178	199	111.9
2	1:2		518	
	1:4	259	254	98.1
	1:8	129	136	105.0
	1:16	65	69	107.4
3	1:2		751	
	1:4	375	373	99.3
	1:8	188	183	97.2
	1:16	94	95	100.8
4	1:2		755	
	1:4	378	397	105.2
	1:8	189	202	107.1
	1:16	94	106	112.1

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (spiked) were stored at -20°C and thawed 5 times, and the mouse TNF- α levels determined. There was no significant loss of mouse TNF- α immunoreactivity detected by freezing and thawing.

13.5.2 Storage Stability

Aliquots of serum samples (spiked) were stored at -20°C, 2-8°C, room temperature (RT) and at 37°C, and the mouse TNF- α level determined after 24 h. There was no significant loss of mouse TNF- α immunoreactivity detected during storage at -20°C and 2-8°C.

13.6 Specificity

The assay detects both natural and recombinant mouse TNF- α . The interference of circulating factors of the immune system was evaluated by spiking these proteins at physiologically relevant concentrations into mouse IL-10 positive serum. There was no crossreactivity detected.

13.7 Expected Values

There were no detectable mouse TNF- α levels found. Elevated mouse TNF- α levels depend on the type of immunological disorder.

14 Ordering Information

+ Europe

Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe
phone: +43 1 796 40 40 ext. 114
fax: +43 1 796 40 40 ext. 400

order@bendermedsystems.com

technical information:

techserv@bendermedsystems.com

+ USA

eBioscience, Inc.
10255 Science Center Drive
San Diego, CA 92121
phone: +1(888) 999 1371
fax: +1(858) 642 2046

contact@ebioscience.com

technical information:

tech@ebioscience.com

www.bendermedsystems.com

www.ebioscience.com

15 Reagent Preparation Summary

15.1 Wash Buffer (1x)

Add **Wash Buffer Concentrate** 20x (50 ml) to 950 ml distilled water.

Number of Strips	Wash Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

15.2 Assay Buffer (1x)

Add **Assay Buffer Concentrate** 20x (5 ml) to 95 ml distilled water.

Number of Strips	Assay Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

15.3 Biotin-Conjugate

Make a 1:100 dilution of **Biotin-Conjugate** in Assay Buffer (1x):

Number of Strips	Biotin-Conjugate (ml)	Assay Buffer (1x) (ml)
1 - 6	0.06	5.94
1 - 12	0.12	11.88

15.4 Streptavidin-HRP

Make a 1:200 dilution of **Streptavidin-HRP** in Assay Buffer (1x):

Number of Strips	Streptavidin-HRP (ml)	Assay Buffer (1x) (ml)
1 - 6	0.03	5.97
1 - 12	0.06	11.94

15.5 Mouse TNF- α Standard

Reconstitute lyophilized **mouse TNF- α standard** with distilled water. (Reconstitution volume is stated on the label of the standard vial.)

16 Test Protocol Summary

1. Determine the number of microwell strips required.
2. Wash microwell strips twice with Wash Buffer.
3. Standard dilution on the microwell plate: Add 100 µl Sample Diluent, in duplicate, to all standard wells. Pipette 100 µl prepared standard into the first wells and create standard dilutions by transferring 100 µl from well to well. Discard 100 µl from the last wells.
Alternatively external standard dilution in tubes (see 9.5.1): Pipette 100 µl of these standard dilutions in the microwell strips.
4. Add 100 µl Sample Diluent, in duplicate, to the blank wells.
5. Add 50 µl Sample Diluent to sample wells.
6. Add 50 µl sample in duplicate, to designated sample wells.
7. Cover microwell strips and incubate 2 hours at room temperature (18° to 25°C).
8. Prepare Biotin-Conjugate.
9. Empty and wash microwell strips 5 times with Wash Buffer.
10. Add 100 µl Biotin-Conjugate to all wells.
11. Cover microwell strips and incubate 1 hour at room temperature (18° to 25°C).
12. Prepare Streptavidin-HRP.
13. Empty and wash microwell strips 5 times with Wash Buffer.
14. Add 100 µl diluted Streptavidin-HRP to all wells.
15. Cover microwell strips and incubate 1 hour at room temperature (18° to 25°C).
16. Empty and wash microwell strips 5 times with Wash Buffer.
17. Add 100 µl of TMB Substrate Solution to all wells.
18. Incubate the microwell strips for about 10 minutes at room temperature (18° to 25°C).
19. Add 100 µl Stop Solution to all wells.
20. Blank microwell reader and measure colour intensity at 450 nm.

Note: If instructions in this protocol have been followed samples have been diluted 1:2 (50 µl sample + 50 µl Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 2).